


PAJUNK®

*The Economic Scissor System
with Disposable Inserts*

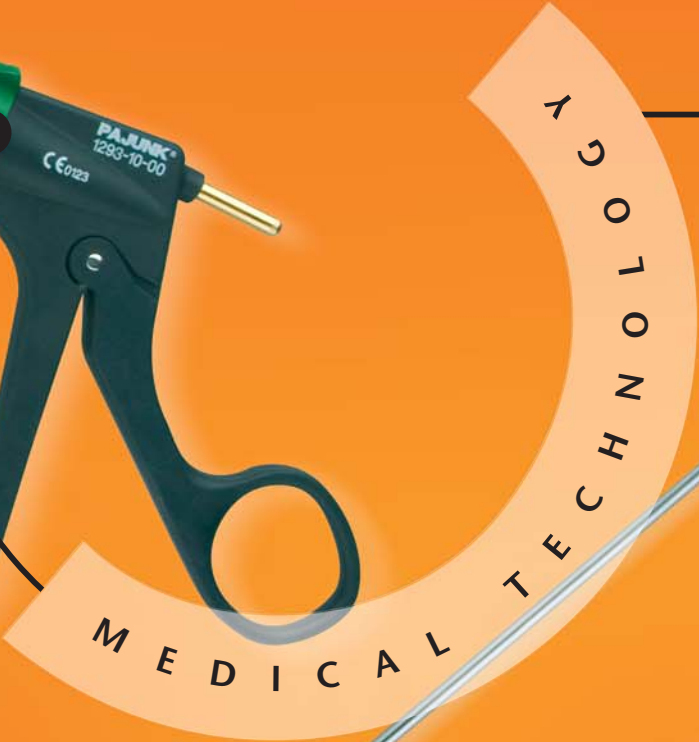

Laparoscopy

Always new, always sharp



The economic scissor system with disposable inserts



PAJUNK® has extended its range of instruments for laparoscopy by a reusable, modular system with disposable inserts, which permits the economic replacement of the sensitive scissor-inserts. And with that, PAJUNK® has created the option to


always work and also to coagulate with new, guaranteed sharp scissors at every intervention.

Three different sizes are available: micro, medium and large. They are packed individually, sterile, and they can be replaced very easily.

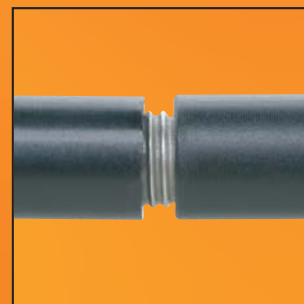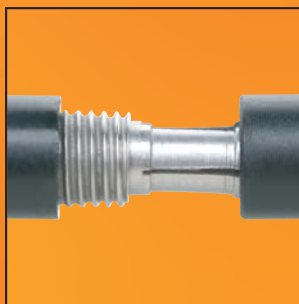

The fast and easy assembly is based on the ingenious locking mechanism. The scissor-insert is clicked into the connecting rod and tightly screwed together with the

guidance tube. Please observe, that the replacement can only be carried out if the scissors are closed.

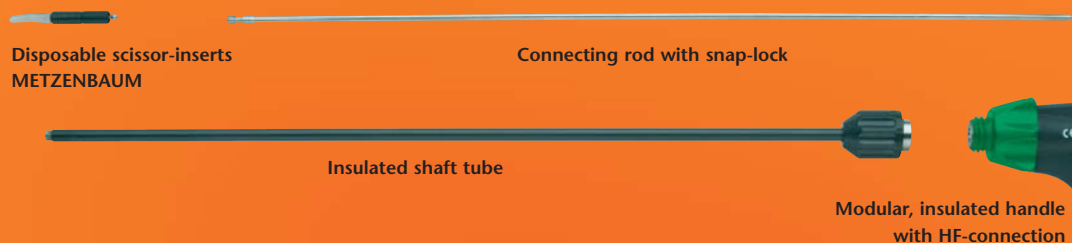

Scissor System with Disposable Inserts

Modular, insulated handle with HF-connection (autoclavable)

	Item no.	PU
	1293-10-00	1




Insulated shaft tube (autoclavable)

Working length 330 mm	1293-15-05	1
Working length 450 mm	1293-16-05	1

Connecting rod with snap-lock (autoclavable)

Working length 330 mm	1293-17-05	1
Working length 450 mm	1293-18-05	1

Disposable scissor – insert METZENBAUM, individually packed, sterile

	Size	Item no.	PU
	large	1293-93-01	10
	medium	1293-93-02	10
	micro	1293-93-03	10

PAJUNK GmbH

Medical Technology

Karl-Hall-Strasse 1

D-78187 Geisingen/Germany

Phone +49 (0) 77 04/92 91-0

Fax +49 (0) 77 04/92 91-6 00

www.pajunk.com