

**PAJUNK®**

*VarioJect Cartridge Syringes  
Made of High-Grade  
Stainless Steel*



Dental



# Dental syringes made of high-grade stainless steel – sturdy, durable and easy to clean

As of now, PAJUNK® VarioJect dental syringes are generally manufactured from stainless steel. The robust material is particularly resistive against damaging and harsh and aggressive cleaning agents. In addition, the special matted surface finish helps to avoid irritating reflections.


The ergonomic design of the syringes permits safe and convenient handling. You can choose from two design variants: with slip-on device or swing-out technology.

When using the slip-on technology, several filled cartridge holders can be prepared and fitted, which guarantees fast exchange during treatment.


In contrast thereto, the swing-out technology permits the cartridge holder to be swung sideways, and only the ampoule itself is replaced.



VarioJect made of stainless steel with swing-out technology

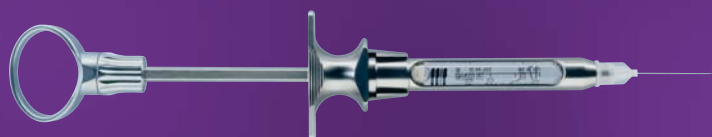

VarioJect made of stainless steel with slip-on device

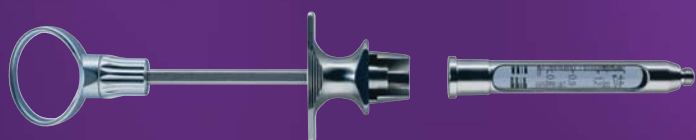

To be used alternatively either with stainless steel cartridge holder, or with single-use, disposable safety cartridge holder VarioSafe



Both techniques support either manual aspiration by means of three extensible aspiration hooklets, as well as self-aspiration. They are suitable for the use of regular commercially available 1.8 ml cylinder ampoules, and for cannulae equipped with a metric thread.

VarioJect made of stainless steel	Item no.
With swing-out technology and cartridge holder for 1.8 ml ampoules, 3 aspiration hooklets and self-aspiration, pcs. per pack: 1 piece	8990-50006
With slip-on device and cartridge holder for 1.8 ml ampoules, 3 aspiration hooklets and self-aspiration, pcs. per pack: 1 piece	8990-60006
Slip-on cartridge holder for 1.8 ml ampoules, with self-aspiration, pcs. per pack: 1 piece	8990-01002
VarioSafe made of stainless steel	Item no.
With slip-on device for VarioSafe single-use, disposable safety cartridge holders with 3 aspiration hooklets, pcs. per pack: 1 piece	8990-41006
VarioSafe single-use, disposable safety cartridge holder with self-aspiration for 1.8 ml ampoules, pcs. per pack: 100 pieces	8990-40002

PAJUNK GmbH Medizintechnologie · Karl-Hall-Strasse 1 · D-78187 Geisingen/Germany  
 Telefon +49 (0) 77 04/92 91-0 Telefax +49 (0) 77 04/92 91-6 00 · www.pajunk.com

Dental